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### User Manual


**Model: HouseSBS-30008**

Thank you for purchasing this House bathroom scale. Please read all operating instructions carefully before use. Scales are precision measuring devices and should always be handled with proper care. To ensure years of reliable service, keep these simple tips in mind:

#### **Battery Installation**

1. This appliance works with battery of 2 pieces AAA ,1.5V.
2. Before using the scale for the first time, check the battery compartment and remove any plastic insulation strips or packaging that may be in place to prevent battery drainage during shipping/storage.
3. If you are replacing the batteries, make sure the positive and negative contacts are properly aligned.
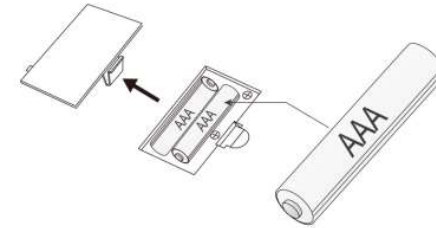
#### **Unit Conversion**

1. The “UNIT” key is on the bottom of the plastic, press the button to select the desired unit, kg=>lb=>lb:st

If the scale is turned off then back on, it will default to the last unit used.

#### **Operation and Weighing**

Step1: Place the scale on a flat, level surface. Do not weigh on carpet as it may affect the accuracy.

Step2: If this is the first time the scale is being used or its position has been

changed, press the center of the the scale to activate the power and wait until the reading returns to zero (0.0) . This is for zero adjustment.

Step3: Step onto the scale platform and stand still. After a few seconds, you will see the weight flashing to indicate a weight reading. Your weight will be displayed for 10 seconds and the scale will automatically turn off.


If you want to measure the weight again, please step off, then step on the platform again.

### **Error Codes**

#### 3. Low Battery Indication:

When the LCD indicate “Lo”, open the battery door and replace the battery with a new one.
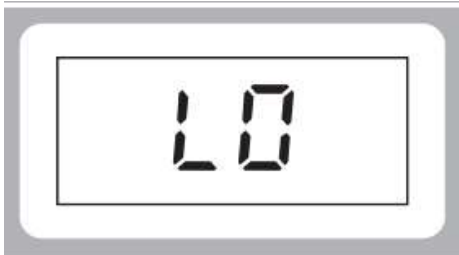
#### 4. Overload Indication:

When the scale is over loaded 180kg, “Err” will appear on the screen. Remove the weight to protect the sensors.

### **Weighing Tips**

- Do not weigh on carpet. Hard floors provide the best weighing surface. Slanted or uneven surfaces should also be avoided.
- For best results, weigh yourself after waking up in the morning. This will give you a more accurate picture of your weight trend. Your weight can fluctuate throughout the day.
- For repeatable results, stand as still as possible when weighing, with your weight centered on the platform.
- Keep the area under and around the scale clear. Inaccuracy can arise from items lodged under the platform.

### **Care Instructions**

- Clean the scale with a slightly damp cloth. Do not use harsh solvents.
- To avoid leakage, remove the batteries when storing the scale for extended periods of time.
- Avoid dropping your scale. The warranty does not cover damage due to rough treatment.
- Check the batteries first if you are having any trouble with your scale. This simple step can remedy most scale issues.
- Store your scale in a clean, dry location. Dust, dirt, and moisture can

accumulate on the weighing sensors and electronics causing inaccuracy of malfunction.

- The scale is for family use only, DO NOT use for commercial purposes.

### Specification

Capacity: 150kg	Graduation: d= 100g
Unit of measure: kg/ lb / lb:st	LCD size: 70x28mm
Auto zero & Auto-off & Auto on	Power: 2x1.5V AAA battery (Not Incl)
Indication: Low battery and overload indication	
Model: House SBS-30008	

This marking indicates that the device must not be disposed of with other household waste. The regulation applies to the entire EU area. To avoid any harm to the environment or health due to improper waste disposal, the device should be properly recycled so that the materials can be reused. If the device will no longer be used, return it to a collection point for waste electronic and electrical equipment or to the shop from which it was purchased. This ensures that the device is recycled in an environmentally safe manner.

### Warranty and liability for defects

“House” brand electrical devices have a two-year warranty from the time of purchase.

The warranty does not cover normal wear and tear of the product or any of its parts, improper handling of the product, such as improper placement and storage, improper connection or installation, damage caused by other external factors, or failure to comply with the operating or maintenance instructions.


Please read the operating instructions carefully as they contain important information.

The product is subject to a statutory liability for defects.

### Note:

1. If the product does not work properly, first check if the fault can be due to other reasons. With electrical devices, the reason is often the interruption of power supply or other incorrect handling.
2. Please provide the following documents and/or information with the defective product:
  - Purchase receipt
  - Product model description/type/brand
  - Detailed description of the defect or problem

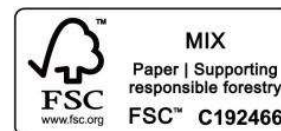In case of a claim for compensation related to the defect, please contact the dealer.

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